







BioLingual® archwires



Size	Profile	Force	Content	 shape 1	 shape 2	 shape 3	 shape 4	 shape 5	 shape 6
.010"	●	20 g	10	270-0125	270-0225	270-0325	270-0425	270-0525	270-0625
.012"	●	30 g	10	270-0130	270-0230	270-0330	270-0430	270-0530	270-0630
.014"	●	40 g	10	270-0135	270-0235	270-0335	270-0435	270-0535	270-0635
.016"	●	50 g	10	270-0140	270-0240	270-0340	270-0440	270-0540	270-0640

BioLingual® archwires – customised packed

Size	Profile	Force	Order No.	Content			
				shape 2	shape 3	shape 4	shape 5
.012"	●	30 g	270-0030	2	3	3	2
.014"	●	40 g	270-0035	2	3	3	2
.016"	●	50 g	270-0040	2	3	3	2



BioCosmetic®-archwires –
durable, perfect aesthetics.

BioCosmetic®-archwires are the first tooth-coloured archwires that retain their white appearance better than powder coated aesthetic archwires and therefore enable long-term aesthetical treatment. They also provide maximum intraoral comfort for the patient, as the white coating covers superelastic NiTi archwires. These Archwires ensure gentle forces and quick treatment due to their extremely low force level.

The intervals for archwire changes with BioCosmetic archwires are the same as for uncoated archwires. They are also suitable for the use of self-ligating brackets, which generally allow longer intervals for archwire changes. The tooth coloured coating of BioCosmetic archwires has a better integrity than powder coating, ensuring perfectly maintained aesthetics. Regardless of whether the bracket slots are made of metal or ceramic, BioCosmetic®-archwires can maintain their white appearance for up to six weeks. In

addition we recommend to use patient friendly brackets with rounded bracket and slot edges as for example the QuickClear® bracket by FORESTADENT. BioCosmetic®-archwires, however, should not be adapted or bent, since sharp objects may damage the white coating.



Description	Order No. Maxillary	Order No. Mandibular	Profile	Force	ø inch	ø mm	Content
BioCosmetic® Superelastic	217-0845	217-0945	●	25 g	.017"	0,43	10
	217-0850	217-0950	●	35 g	.019"	0,48	10
	217-1440	217-1540	■	55 g	.016" x .016"	0,40 x 0,40	10
	217-2040	217-2140	■	90 g	.016" x .022"	0,41 x 0,56	10
	217-2046	217-2146	■	135 g	.018" x .025"	0,46 x 0,64	10

FORESTADENT Bernhard Förster GmbH
Westliche Karl-Friedrich-Str. 151 • 75172 Pforzheim (Germany)
Phone: +49 (0) 7231 459-0 • Fax: +49 (0) 7231 459-102 • info@forestadent.com

www.forestadent.com

FORESTADENT®
GERMAN PRECISION IN ORTHODONTICS

The heart piece of your treatment
Our Forestadent BioArchwires

www.forestadent.com

FORESTADENT®
GERMAN PRECISION IN ORTHODONTICS

The intelligent archwire for every treatment stage



BioStarter®: Round archwires for a perfect start to your orthodontic treatment

BioTorque®: Square and Rectangular archwires for efficiency

when you need to express detail

BioLingual®: Round archwire for intelligent lingual treatment





BioCosmetic®: Durable, perfect aesthetics

The correct force for each treatment phase.

In orthodontics, just like in car racing, the correct application of force and torque is crucial: too much of either and you will get unwanted side effects, too little and nothing happens at all. Archwires are the engine of orthodontic tooth movement. During each treatment phase the correct choice of archwire should guarantee the correct force application. You will reach your treatment goals faster without unwanted side effects. When using BioArchwires from FORESTADENT you can choose the correct force levels for any stage of your treatment.

Faster treatment progress possible due to reduced friction.

The uniquely smooth surface of BioArchwires reduces the friction between archwire and bracket by up to 30 %. Due to the combination of the unique characteristics these intelligent wires from FORESTADENT work gentle and efficient. The orthodontic treatment can be completed faster. When using BioArchwires treatment can potentially be completed with fewer archwires in a shorter period of time – a more economical solution for the clinician and the patient.

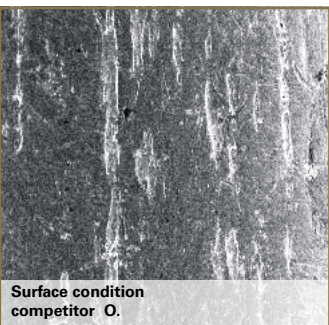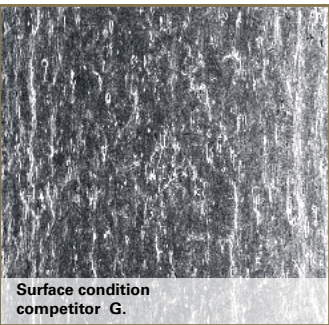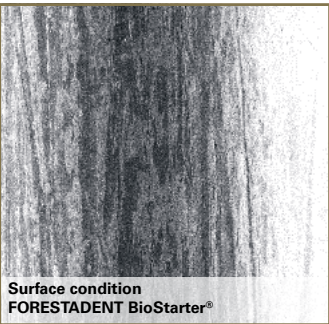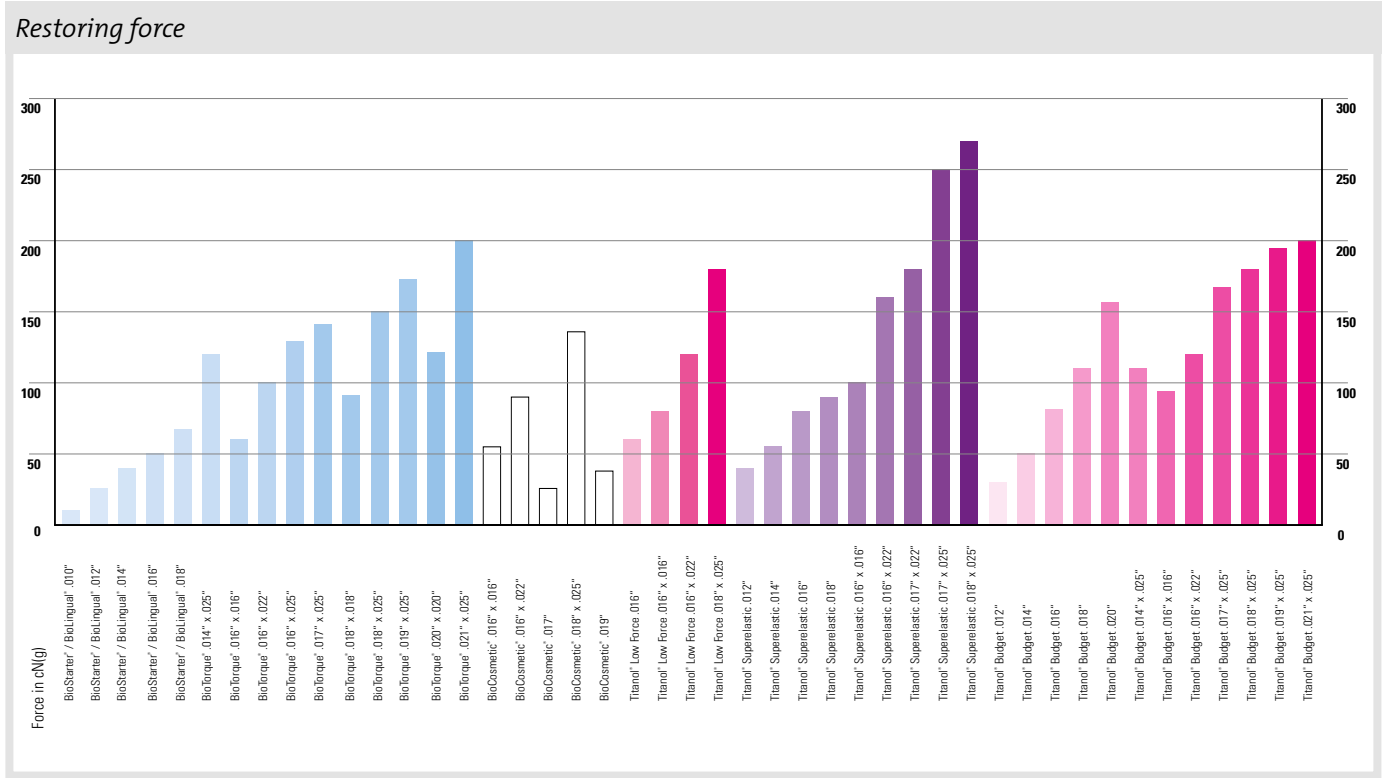

The uniquely smooth surface of BioStarter® and BioTorque® archwires reduces friction and allows faster and more efficient treatment progress.

BioStarter® round 0.010" – for particularly sensitive patients

The smallest BioStarter® archwire with the dimension 0.010" broadens the application of BioStarter® archwires even further. This archwire will work in the most adverse conditions, e.g. periodontally compromised patients, and malocclusions can be corrected with the lightest of forces and extremely gentle. Until now stranded wires with high friction and poor elasticity had to be used for these applications leading to extended treatment times. The BioStarter® archwire is available in force levels also suitable for sensitive and anxious patients.

BioTorque® square and rectangular wires – efficiency in every direction

Torquing of teeth with rectangular archwires generates intensive bone remodeling and high compressive forces that can damage the dentition. The thermoactive BioTorque® archwires from FORESTADENT can reduce the risk of root resorption due their low forces. Due to the extremely precise dimensions of the FORESTADENT wires exact torquing is a lot easier. BioTorque® archwires feature the same unique surface finish as the BioStarter® archwires helping to keep friction to a minimum. Together with the extreme superelasticity BioTorque® archwires work very efficiently and super fast. BioTorque® archwires can be used in difficult cases, e.g. severe crowdings or rotations.

BioLingual® round archwire for intelligent lingual treatment

Thermoactive bioarchwires are not only available for vestibular but also for lingual treatment. The BioLingual® archwires are mushroom shaped and can be ordered in six different shapes with four dimensions each (.010", .012", .014", und .016"). The shape, optimized for lingual application, was developed in cooperation with orthodontic experts. The nickel-titanium BioLingual® archwires generate their optimal force at body temperature, in the same way as the popular BioStarter® and BioTorque® archwires.

What are you waiting for?

It's time for a change. Order your BioStarter®, BioTorque® and BioLingual® wires today and change to the intelligent and patientfriendly technology from FORESTADENT. FORESTADENT technology – for your and your patients benefit.




River Finish Gold

Over 111 years of working with surface polishing is ideal for creating the optimum arch wires. Inspired by the smooth surface of river stones, FORESTADENT archwires undergo a unique treatment process Called River Finish Gold. The result is ultra-smooth archwires with low friction. The best possible foundation for a perfect treatment result.

BioStarter®, BioTorque®

Individually packed

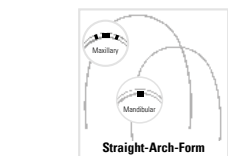


Straight-Arch-Form



Euro-Smile-Form

Description	Order No. Maxillary	Order No. Mandibular	Order No. Maxillary	Order No. Mandibular	Profile	Force	ø inch	Content
BioStarter®	203-0825	203-0925	203-1825	203-1925	●	20 g	.010"	10
	203-0830	203-0930	203-1830	203-1930	●	30 g	.012"	10
	203-0835	203-0935	203-1835	203-1935	●	40 g	.014"	10
	203-0840	203-0940	203-1840	203-1940	●	50 g	.016"	10
	203-0845	203-0945	203-1845	203-1945	●	70 g	.018"	10

BioTorque®	203-2035	203-2135	203-4035	203-4135	■	120 g	.014"x.025"	10
	203-1440	203-1540	203-3440	203-3540	■	60 g	.016"x.016"	10
	203-2040	203-2140	203-4040	203-4140	■	100 g	.016"x.022"	10
	203-2440	203-2540	203-4440	203-4540	■	135 g	.016"x.025"	10
	203-2044	203-2144	203-4044	203-4144	■	140 g	.017"x.025"	10
	203-1445	203-1545	203-3445	203-3545	■	90 g	.018"x.018"	10
	203-2046	203-2146	203-4046	203-4146	■	150 g	.018"x.025"	10
	203-2048	203-2148	203-4048	203-4148	■	170 g	.019"x.025"	10
	203-1450	203-1550	203-3450	203-3550	■	130 g	.020"x.020"	10
	203-2053	203-2153	203-4053	203-4153	■	200 g	.021" x .025"	10

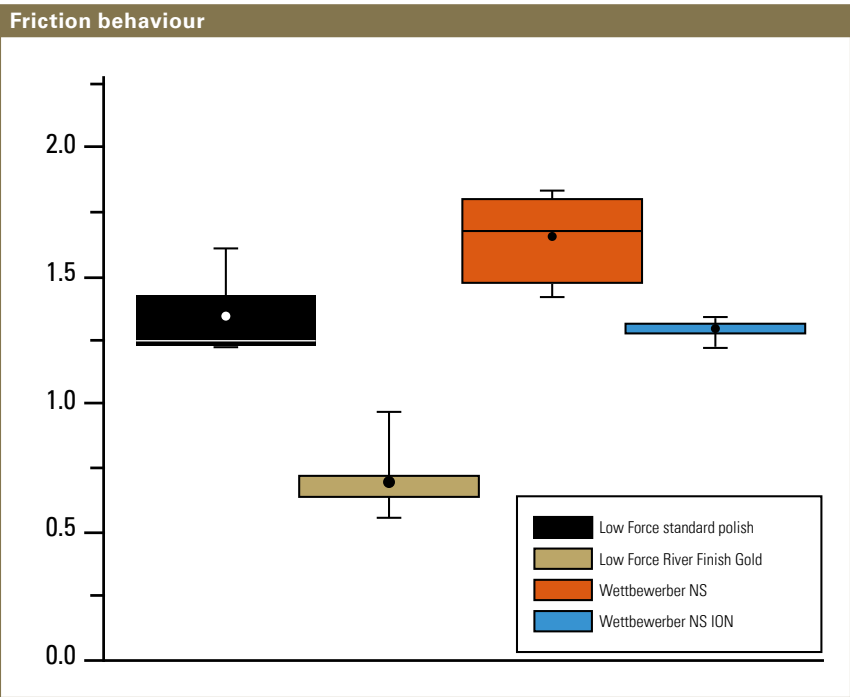

The effect of different surface coatings (River Finish Gold and ion implantation) on the friction behaviour in comparison with uncoated archwires.